


PULSED
LIGHT



Unforgiving with hair and gentle on the skin

With Modula hair removal becomes more effective and painless: the dual cooling technology maximizes product performance, and of course the comfort of the client. The efficiency of the water system replaces the obsolete air system, thus allowing the achievement of higher energy values (35 J/cm²) and acting more decisively on hair or skin spots. The SCS Skin Cooling System operates on contact cooling the skin and preventing discomfort and redness. Client safety always comes first. Wavemed designed the **FCC Filter Check Control**, a device that terminates the treatment if the filter is not inserted thereby avoiding accidents and damage to the skin.

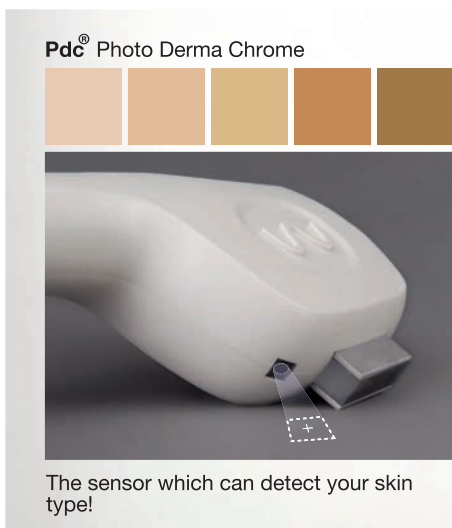

Treatments tailored for each skin type

With a single gesture, you can personalize the treatment and precisely calibrate the gauges for every client. Thanks to the innovative PDC Photo Derma Chrome technology, a photo sensor on the hand-piece records the exact colour shade of the skin and the system automatically selects the most suitable programme for that skin type.



A variety of six treatments in a single hand-piece

Following the logic of practicality and minimum size, Modula uses interchangeable filters. This means that a single hand-piece is able to perform six different treatments. By selectively cutting the broad spectrum light beam, we can obtain the selected wavelengths for specific actions: from hair removal to photo-rejuvenation, from acne therapy to vascular treatments.

XL coverage

To shorten the work time, the hand-piece comes in a version with an extra-large coverage that allows you to work on a large area (6 cm²), decreasing the duration of each session, thus lowering costs and increasing the pleasure.

TREATMENT AREA
50 x 12 mm

The results can be seen Speed as an asset.

With its multi-action mode and minimum recharge times between pulses (about 1 second), Modula bombards the hair follicles in a much more intense and effective manner. Permanent hair removal is no longer a mirage and becomes an achievable goal in a few sessions.

It all fits on a trolley

The well designed details reduce the size to a minimum, so that if necessary you can store it in its trolley and it can be transported easily, according to your needs. This fast, powerful device can be transported as hand baggage, able to offer to you quality and convenience anywhere.

Details that make you feel safe

Wavemed has thought of everything, even the smallest details. A complete system of sensors controls all aspects of technology and notifies you whenever necessary, from the level of lamp life remaining to the level of the flow of the hydraulic system or to that of the water in the cooling system.

Modula

PULSED LIGHT

General	FCC Filter Check Control PDC Photo Derma Chrome Display Touch Screen 10,4"
---------	--

Power source	110 – 240 VAC 50-60 Hz
Absorption	1200 VA
Device type	BF
Medical device	Class IIb - CE 1936
Casing	Fire proof ABS
Insulation class	I
Fuses	2 x T6,3A
Fluency	26 J/cm ² (aesthetics)
Fluency	35 J/cm ² (medical)
Area spot	50 x12 mm
Lamp shots	100 000
Charge time	1 sec
Burst mode	5 consecutive pulses
Interchangeable filters	640-1200 nm hair removal 600-1200 nm hair removal (optional) 570-1200 nm photo-rejuvenation 550-1200 nm deep vascular treatment (optional) 520-1200 nm surface vascular treatment (optional) 400-1200 nm acne
Cooling	Purified water
Skin cooling	Cooled crystal
PDC - Photo Dermal Control	Automatic skin type recognition
Connectivity	Wifi , LAN, USB
Size l x h x d	400 x 310 x 230 mm
Weight	20 kg
Accessories	Cart, trolley, additional filters
Accuracy of displayed values	± 20 %

WAVEMED SRL, Via Campo Gillaro, 29A, 00030 San Cesareo (RM) - ITALY
Phone +39 06 97652423 - Fax +39 06 97652424
www.wavemed.it - info@wavemed.it



Management System
ISO 9001:2008
EN ISO 13485:2012

www.tuv.com
ID 9105085995